

Product Profile

Laboratory instruments for quality control, analysis and calibration



Stanhope-Seta specialise in the design and manufacture of laboratory test instruments which are widely used to measure product quality and consistency. From its manufacturing base in the United Kingdom the company exports high quality instrumentation worldwide.

The wide range of instruments, manufactured in the UK are regarded as a worldwide benchmark for quality control, with products that offer advanced performance, accuracy and reliability. The company is certified to ISO 9001 and ISO 14001. Stanhope-Seta is also proud to hold a King's Award for Enterprise in the International Trade category and four Queen's Awards for Industry in the Innovation and International Trade categories.



In addition, Seta produces a comprehensive range of Reference Materials including Certified Reference Materials (CRMs), Multi Test Verification Materials (MTVMs) and Single Test Verification Materials (STVMs), including those produced under ISO 17025 and ISO 17034 accreditation.



Stanhope-Seta is an active member of major standardisation bodies including ASTM, IP, ISO, BSI and DIN. By working closely with industry our research and design group actively contributes to test method development and product specifications.



Contents

Air Release..... 22

Biodiesel and Ethanol Content 14-15

Carbon Residue and Total Sediment 20

Cetane Number 30

Cloud and Pour Point 21

Crude Oil Testing, and Sediment by Filtration 26-27

Density 13

Dielectric Measurement 28

FAME Measurement 29

Filter Blocking Testing 11

Flash Point..... 6-10

Grease Testing..... 24

Gum Residue, Particulates by Filtration 16

H₂S Analysis..... 31

Oxidation, Silver and Copper Corrosion 17-19

Particle Counting..... 32-36

Reference Materials..... 38-39

Vapour Pressure 12

Viscosity..... 25

Water Separability..... 23

Water Separation, WSI 37



Application	Jet A-1 DefStan 91-091	Aviation Turbine Fuel ASTM D1655	Aviation Turbine (Synthesized) D7566	Aviation Gasoline ASTM D910	Motor Gasoline EN 228	Motor Gasoline ASTM D4814	Diesel EN 590	Diesel ASTM D975
Ash							EN ISO 6245	ASTM D482
Carbon Residue							EN ISO 10370	
Cetane Number							IP 617; EN 17155	ASTM D8183
CFPP							EN 116	ASTM D6371
Cloud Point							EN 23015	ASTM D2500
Colour	ASTM D156; D6045	ASTM D156; D6045	ASTM D156; D6045					
Copper Corrosion	IP 154; ASTM D130	IP 154; ASTM D130	ASTM D130; IP 154	ASTM D130	EN ISO 2160	ASTM D130	EN ISO 2160	ASTM D130
Silver Corrosion		IP 227				ASTM D7667; D7671		
Dielectric								
Distillation	ASTM D86; IP 123	ASTM D86	ASTM D86	ASTM D86	EN ISO 3405	ASTM D86	EN ISO 3405	ASTM D86
Existent Gum	ASTM D381; IP 540	IP 540; ASTM D381	IP 540; ASTM D381		EN ISO 6246	ASTM D381		
FAME Content	ASTM D7797; IP 583	ASTM D7797; IP 583					ASTM D7963; ASTM D8274; IP 639	ASTM D7963; ASTM D8274; IP 639
Filter Blocking Tendency							IP 387; 618	
Freeze Point	ASTM D2386; IP 16	ASTM D2386	D2386	ASTM D2386				
Flash Point	ASTM D56; D93; D3828; D7236; IP 34; IP 170; IP 523; IP 534	ASTM D56; D93; D3828; D7236; IP 170; IP 523; IP 534	ASTM D56; D3828; D7236; IP 170; IP 523; IP 534				EN ISO 2719	ASTM D93; D3828
Hydrogen Sulfide								
Oxidation Stability				ASTM D873	EN ISO 7536	ASTM D525	EN ISO 12205	ASTM D2274
Particle Counting	ASTM D7619; IP 565	IP 565; ASTM D7619					IP 630	ASTM D7619
Particulates Contamination (Filtration)	IP 423; ASTM D5452	ASTM D2276; D5452; IP 216; 423	ASTM D2276; D5452; IP 216; 423				EN 12662	ASTM D6217
Pour Point								
Total Sediment								
Unconverted Fatty Acids	ASTM D7797; IP 583							
Vapour Pressure		ASTM D323; D4953; D5190; D5191; IP 69; 394		ASTM D323; D5191	EN 13016-1	ASTM D5191; D4953; D5482; D5188; D6378		
Viscosity	ASTM D445; IP 71	ASTM D445; IP 71	ASTM D445; IP 71				EN ISO 3104	ASTM D445
Water Separation	ASTM D8073; IP 624	ASTM D8073; IP 624						

Laboratory instruments for quality control, analysis and calibration



	Marine Fuel ISO 8217	Fuel Oil ASTM D396	Bio-Heating Fuels EN 14213	FAME EN 14214	Bio-Fuel Blend Stock B20 ASTM D7467	Bio-Fuel Blend Stock B100 ASTM D6751	Instrument	Seta Reference Number
	ISO 6245	ASTM D482	ISO 3987	ISO 3987	ASTM D482	ASTM D874	Ash Furnace	99220-2
	ISO 10370; IP 398		EN ISO 10370	EN ISO 10370		ASTM D4530; D189	Micro Carbon Residue Tester	97400-3
				IP 617; EN 17155	ASTM D8183	ASTM D8183	AFIDA Cetane Number Analysis	SA6000-0
	IP 309		EN 116	EN 116	ASTM D6371		Cloud & Pour Point Bath	94160-0
	ISO 3015	ASTM D2500		EN ISO 23015	ASTM D2500	ASTM D2500	Cloud & Pour Point Bath Cloud & Pour Point Cryostat	94160-0 93531-8
							Colorimeter	15260-4 15320-2
		ASTM D130		ISO 2160; IP 154	ASTM D130	ASTM D130	Copper and Silver Corrosion Baths	11310-0 11405-0
							Copper and Silver Corrosion Baths	11310-0 11405-0
							JetDC	88500-0
		ASTM D86			ASTM D86		Setastill Distillation	11860-3
							Existent Gum Bath	12200-3 12210-0
	ASTM D7963; ASTM D8274; IP 639	ASTM D7963; ASTM D8274; IP 639	ASTM D8274; IP 639		ASTM D8274; IP 639		SetaCheck BioDiesel FIJI Instruments	SA5550-0 SetaCheck SA5000-2 - FIJI Jet SA5100-0 - FIJI Distillates SA5200-0 - FIJI Multi Fuel
							Multi-Filtration Tester Cold Filter Blocking Tester	91600-4 91670-2
							Freezing Point Apparatus	16990-2
	EN ISO 2719	ASTM D93; D3828	EN ISO 3679; IP 523; 524	EN ISO 3679; 2719; IP 523; 524	ASTM D93; D3828; D56	ASTM D3828; D93	Pensky-Martens Setaflash Series 3 Range Setaflash Series 8 Range	35000-0 Series 3 Series 8
	IP 570						H ₂ S Analyser with VPP	SA4000-4
	ISO 12205						Oxidation Bath TOST Bath TOST Dry Block Bath	16900-7 16640-2 16645-2
					ASTM D7619		AvCount3 AvCount Lite	SA1100-0 SA1800-2
			EN 12662	EN 12662	ASTM D6217	ASTM D6217	Filtration Tests	Contact Seta for Filtration
	ISO 3016	ASTM D97	ISO 3016				Cloud & Pour Point Bath Cloud & Pour Point Cryostat	94160-0 93531-8
	ISO 10307						Total Sediment Tester	16120-2
							FIJI Instruments	SA5000-2 - FIJI Jet SA5200-0 - FIJI Multi Fuel
							SetaVap2 SetaVap4	81000-2 80600-0
	EN ISO 3104; IP 71	ASTM D445	EN ISO 3104; EN ISO 3105; IP 71	EN ISO 3104	ASTM D445	ASTM D445	KV-6 Viscometer Bath KV-2 Viscometer Bath	84200-3 94710-4
							WSI Analyser	SA9000-0

Setaflash® Series 3 Small Scale Flash Point Testing

The safest choice

Approved in over 1000 international product specifications and regulations


ASTM D3278; ASTM D3828; ASTM D4206; ASTM D7236; ASTM D8174; E502; IP 523; IP 534; IP 602; ISO 3679; ISO 3680 (obs); ISO 9038; EPA 1020 B; CLP Regulations

The versatile range of Setaflash® Series 3 instruments are ideal for use in the laboratory, production line or for portable test applications.

- Electric ignitor or gas (model dependant)
- Test time of under 2 minutes
- Small sample size, 2 or 4 mL
- Temperature range 10 to 135 °C or ambient to 300 °C
- Portable, lightweight, compact design
- Suitable for unknown samples using ramp mode
- Automatic flash detection
- Automatic barometric correction
- Full touch screen display
- USB output and results storage
- Sub-Ambient tests and rapid cool-down



Small Scale CRM Standards:	Nominal Value:	Volume:	ISO 17034
99878-3	75 °C	50 mL	✓
99879-0	192 °C	50 mL	

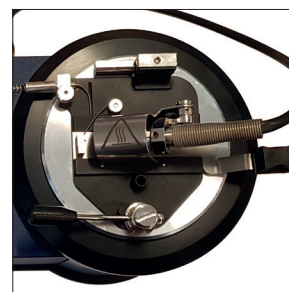

	Series 3 Closed Cup 30000-3	Series 3e Closed Cup 30020-0	Series 3 Open Cup 31000-0	Series 3 ActiveCool 33200-3	Series 3e ActiveCool 33220-0	Series 3 ActiveCool (corrosion resisting) 33250-3	Series 3e ActiveCool (corrosion resisting) 33270-0
Temperature range	Ambient to 300 °C	Ambient to 300 °C	Ambient to 300 °C	10 to 135 °C	10 to 135 °C	10 to 135 °C	10 to 135 °C
Ignitor	Gas	Electric	Gas	Gas	Electric	Gas	Electric
Cup material	Aluminium	Aluminium	Aluminium	Aluminium	Aluminium	Stainless steel insert	Stainless steel insert
Ramp rate	up to 6 °C/min	up to 6 °C/min	n/a	up to 6 °C/min	up to 6 °C/min	up to 6 °C/min	up to 6 °C/min
Heating/cooling method	Cartridge	Cartridge	Cartridge	Peltier cell	Peltier cell	Peltier cell	Peltier cell
Results download	USB	USB	n/a	USB	USB	USB	USB
Size (HxWxD) / Weight	19.5x29.5x14 cm / 3 kg	19.5x29.5x14 cm / 3 kg	19.5x29.5x14 cm / 3 kg	19.5x29.5x14 cm / 3 kg	19.5x29.5x14 cm / 3 kg	19.5x29.5x14 cm / 3 kg	19.5x29.5x14 cm / 3 kg

Setaflash® Series 8 Small Scale Flash Point Testing

The safest choice

Approved in over 1000 international product specifications and regulations

ASTM D3278; ASTM D3828; ASTM D7236; ASTM D8174; E502; IP 523; IP 534; ISO 3679; ISO 3680 (obs);
EPA 1020 B; CLP Regulations

The Setaflash® Series 8 range are automated closed cup flash point testers with enhanced functionality over a wide temperature range.

- Automatic dipping and flash detection
- Electric ignitor (with gas option)
- Test time of under 2 minutes
- Small sample size, 2 or 4 mL
- Temperature range -30 to 135 °C or ambient to 300 °C
- ActiveCool electronic Peltier cooling, no external water bath
- Suitable for unknown samples using ramp mode
- Automatic barometric pressure correction
- Full touch screen display
- Fire detection (relay and alarm)



Small Scale CRM Standards:	Nominal Value:	Volume:	ISO 17034
99878-3	75 °C	50 mL	✓
99879-0	192 °C	50 mL	



	Series 8 High Temperature 82000-2	Series 8 High Temperature 82050-2	Series 8 ActiveCool 82100-2	Series 8 ActiveCool 82110-2	Series 8 ActiveCool (corrosion resisting) 82150-2	Series 8 ActiveCool (corrosion resisting) 82160-2
Temperature range	Ambient to 300 °C	Ambient to 300 °C	-30 to 135 °C	-30 to 135 °C	-30 to 135 °C	-30 to 135 °C
Ignitor	Electric	Gas	Electric	Gas	Electric	Gas
Cup material	Aluminium	Aluminium	Aluminium	Aluminium	Stainless steel insert	Stainless steel insert
Ramp rate	up to 10 °C/min	up to 10 °C/min	up to 10 °C/min	up to 10 °C/min	up to 10 °C/min	up to 10 °C/min
Heating/cooling method	Ceramic hot plate and fan	Ceramic hot plate and fan	Peltier cell	Peltier cell	Peltier cell	Peltier cell
Results download	USB/RJ45	USB/RJ45	USB/RJ45	USB/RJ45	USB/RJ45	USB/RJ45
Size (HxWxD) / Weight	30x34x38 cm / 8 kg	30x34x38 cm / 8 kg	30x34x38 cm / 8 kg	30x34x38 cm / 8 kg	30x34x38 cm / 8 kg	30x34x38 cm / 8 kg

SETAFLASH CCC (35600-0)

ASTM D7094; ASTM D6450; IP 620





The Setaflash CCC is a state of the art Continuously Closed Cup flash point tester designed to deliver fast, accurate and reliable results with minimal sample volume. It combines advanced safety features, automation and a user-friendly touchscreen interface to streamline flash point testing. The Setaflash CCC offers unmatched performance in a compact, modern design.

- 2 mL sample size, reducing waste and hazard
- Electric spark ignition
- Real time display of test progress
- Automatic flash detection
- Removable cup, easy to clean
- Automatic magnet stirrer optimises operation and test validity
- Custom methods and profiles, including search
- Large, intuitive full colour touchscreen display
- Automatic barometric correction
- Internal results storage
- LIMS, network and USB connectivity
- Statistical quality control software (SQC)



Temperature range	-5 °C (-15 °C with external chiller) to 420 °C
Ignition system	Electric
Flash detection	Pressure increase >20kPa
Ramp rate	Standard rates - 2.5°C (D7094) , 5.5°C (D6450)
Fire detection	Thermal
Sample size	1 mL (ASTM D6450), 2 mL (ASTM D7094)
Cooling method	Peltier and Liquid cooled
Result downloads	LIMS, USB, Printer
Size (HxWxD) / Weight	430x 249x320 / 16.3 KG

To further enhance laboratory efficiency, the Setaflash CCC features a quick test function that pre-heats the instrument to the desired temperature before the sample is introduced, minimising downtime and maximising throughput. After testing, the automatic cooling mode rapidly reduces the temperature using Peltier assisted cooling, ensuring the instrument is ready for the next sample in the shortest possible time.

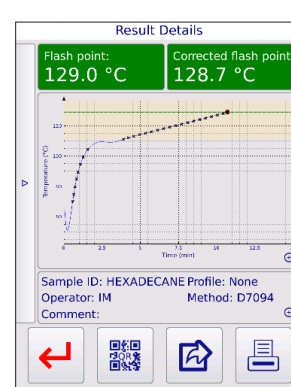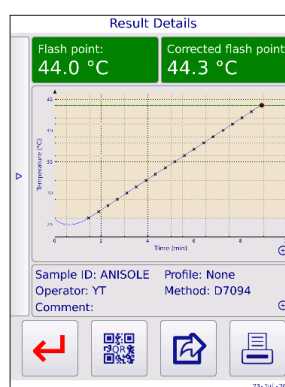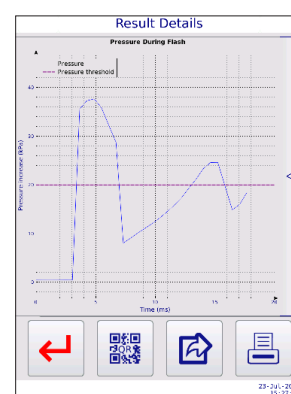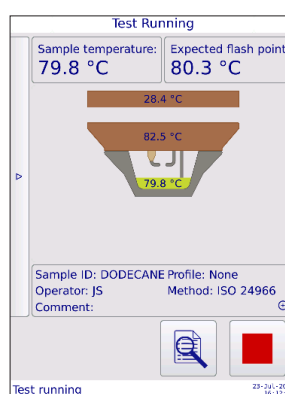

PM-93 Pensky-Martens Closed Cup Flash Point Tester (35000-0)





ASTM D93; D8175; IP 34; ISO 2719 Procedures A, B and C

The Seta PM-93 provides operators with high level functionality combined with class leading ease of use, robustness and safety.

- Fast, accurate and safe
- Easy operation
- Single action lifting pod
- Unique SafeFlash fire extinguishing system
- Seta Ignite - a robust and long lasting ignitor
- 30 programmable test profiles, test methods and sample names
- Large touch screen
- Memory storage for 2000 results
- Statistical quality control software (SQC)



Cooling	Forced air (integral fan)
Ignition system	Electric hot-wire or gas flame
Flash detection	Thermocouple
Heating rate	Fast heating mode (>15 °C/min) and standard rate; 5.5 °C/min; 3 °C/min; 1.3 °C/min; 1 °C/min
Temperature range	Ambient +5 to 400 °C (Ambient +41 to 752 °F)
Barometric pressure correction	Automatic correction with built-in pressure sensor
Fire detection	Thermal. Optional integral inert gas Fire Extinguisher (35002-0)
Size (HxWxD) / Weight	39x24x50 cm / 25 kg

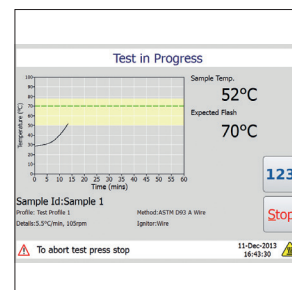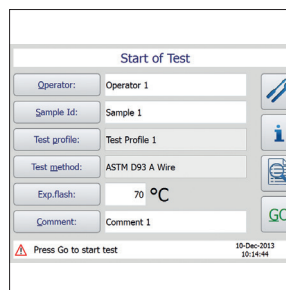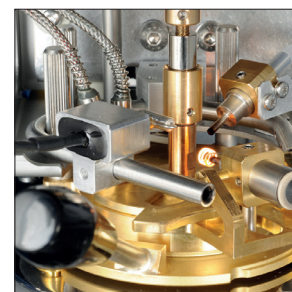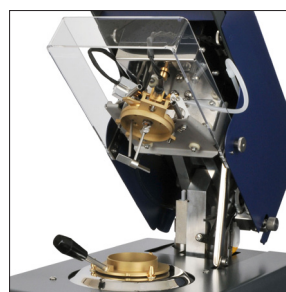

The PM-93 incorporates a heated air bath, a removable sample cup with stirrer, a removable lid with automatic shutter, interchangeable gas and electric hot wire ignitors and a flash detector. The test process is fully automated. Once a sample is added, test parameters are entered and the test is initiated, the test proceeds without any further operator intervention.

At the end of the test, the PM-93 displays the test result, which it automatically corrects for variations in atmospheric pressure using the integral barometer. Test results are saved in the system memory and can be automatically printed, emailed and exported to a USB device or LIMS.

The built-in cooling system helps maximise test efficiency, cooling the sample cup rapidly so that another test can start as soon as possible.




Pensky-Martens CRM Standards:	Nominal Value:	Volume:	ISO 17034
SETA-0412-0051	75 °C	100 mL Pack of 3	✓
SETA-0412-0052	140 °C	100 mL Pack of 3	✓
SETA-0412-0053	225 °C	100 mL Pack of 3	✓

Seta C-92 (35300-0)

ASTM D92; ASTM D8254; IP 36; IP 403 (obs); ISO 2592; DIN 51 376; NF T60-118; JIS K 2265-4; AASHTO T48

The Seta C-92 is an automated Cleveland open cup flash and fire point tester. Featuring the most advanced, built in modern software it is both flexible and easy to use.

- Easy operation
- Electric hot-wire or gas flame
- Real time display of test progress
- Custom methods and profiles
- Full touch screen display
- Automatic correction for atmospheric pressure
- Rapid cooling system
- Automatic snuffer
- Internal result storage
- LIMS, network and USB connectivity
- Statistical quality control software (SQC)
- Single action draught screen

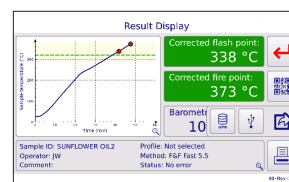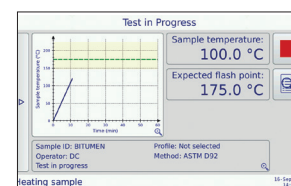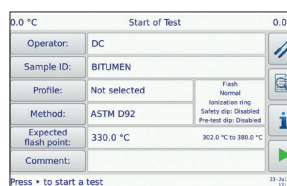

Temperature range	Ambient to 400 °C
Ignition system	Electric hot-wire or gas flame (field switchable)
Flash detection	Ionisation ring (double)
Ramp rate	Automatic or custom settings up to 17 °C/min
Fire detection	Inert gas fire extinguisher and detector
Cooling method	Rapid forced air (integral fan)
Results download	Memory storage for 500,000 test results LIMS, Ethernet, USB
Size (HxWxD) / Weight	26x37x45 cm / 20 kg

Cleveland CRM Standards:	Nominal Value:	Volume:	ISO 17034
SETA-0412-0066	150 °C	100 mL Pack of 3	✓
SETA-0412-0067	265 °C	100 mL Pack of 3	✓
SETA-0412-0068	350 °C	100 mL Pack of 3	✓



Seta Multi Filtration Tester (MFT) (91600-4)


ASTM D2068; ASTM D8386; IP 387; IP PM EA/13; CGSB-3.0 No.142

The Seta MFT is a fully automated instrument, designed to test the Filter Blocking Tendency (FBT) of diesel, biodiesel (B100 & B5/7/20/30), gas oil, gas turbine fuel and kerosene.

- Fully automated test
- Interlaboratory study proven with ASTM and EI validated precision
- Large colour touch screen display
- Graphical display of sample flowrate
- Results in less than 15 minutes
- Automatic calculation of FBT and EFBT
- Statistical Quality Control function
- Storage of test results
- Ethernet and USB port
- RS232 interface for printing
- Simple calibration/verification



Filter blocking tendency range	1.0 to 30 (low number is best)
Sample size	300 mL
Maximum pressure	200 kPa
Flow rate	20 ml/min, adjustable
Power	110/240 V, 50/60 Hz, 300 W
Size (HxWxD) / Weight	35x30x40 cm / 9.8 kg

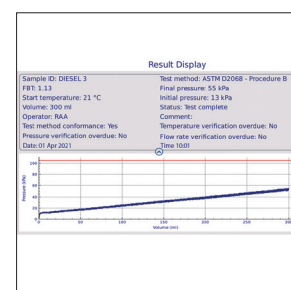

FBT Standards:	Nominal Value:	Volume:
91668-0	1.5 to 3.5 FBT	350 mL

Seta Cold Filter Blocking Tester (91670-2)

IP 618; BS EN 590 National Annex; IP 387; ASTM D2068 Procedure B

The Seta CFBT is a benchtop unit which measures the cold filter blocking tendency (CFBT) and cold filterability of middle distillate fuels containing biodiesel.

CFBT Range	1.0 to 30
Flow rate	20 ml/min, adjustable in software
Sample size	750 mL
Temperature range	-5 °C to +40 °C
Size (HxWxD) / Weight	32x53x54 cm / 24 kg



SetaVap4 Automatic Vapour Pressure Tester (80600-2)

ASTM D5191; ASTM D5188; ASTM D6377; ASTM D6378; ASTM D8009; EN 13016-1 & 3; IP 394
Correlates with ASTM D323; ASTM D1267 ASTM D2533; ASTM D4953; ASTM D5482

Gasoline • Crude Oil • LPG

The SetaVap4 provides users with the latest generation of automated moveable piston based technology for vapour pressure testing.

- Piston design for triple and single expansion test methods
- Peltier cells for heating and cooling
- Integrated shaker for crude oil measurements with variable frequency control
- Large LCD touchscreen display
- Simple test procedure
- Internal results storage
- Compact, rugged design ideal for portable field or laboratory use
- DVPE, RVP, EPA and CARB calculations
- Statistical quality control software
- Optional sample temperature probe
- Optional SETA Digital Density with DensAI Technology, ASTM D4052



Pressure range	0-1000 kPa
Pressure resolution	0.001 kPa
Temperature range	0 to 120 °C (no external cooling required)
Temperature stability	±0.01 °C
Sample volume	1 mL, plus rinsing
Measurement time	approx 5 mins for D5191
Vapour to liquid ratio	0.02:1 to 20:1
Sample introduction	Automated via built in piston, Transfer tube or syringe
Size (HxWxD) / Weight	39x24x30 cm / 9.5 kg



Shown with optional temperature probe

Vapour Pressure Standards:	Test Method:	Range:	Volume:
80610-0	ASTM D5191; EN 13016-1; ASTM D6378; EN 13016-3	P _{TOT} 111.6 to 114 kPa VP ₄ 106.7 to 109.1 kPa	200 mL

Crude Oil Manual Piston Cylinder (80615-0)

The Manual Piston Cylinder (MPC) enables crude oil samples to be taken from a storage tank, rail car or pipeline without loss of volatile light ends for a safe and accurate test.

- Up to 1380 kPa (200 psi)
- 500 ml sample volume
- Dual scale gauge (psi/kPa)



SetaVap4 Vapour Pressure with Digital Density (80607-2)

ASTM D4052; IP 365; ISO 12185

Seta Digital Density is an additional test capability for SetaVap4, providing combined density and vapour pressure testing in a single unit (80607-2).

- Accurate and precise density measurement
- Fully compliant to ASTM D4052
- Automatic Bubble Detection using our patented DensAI technology
- Density measurement via glass oscillating U-tube
- User friendly density calibration screen

The density sensor uses a temperature and humidity controlled oscillating U-tube, allowing for consistent performance across all laboratory testing conditions. The combination of our patented technology, DensAI and a built-in HD camera allows for seamless automatic bubble detection.

Date/Time: 18 Mar 2025 11:45	Comment: BATCH A
User: RAA	Cell temperature: 20.0 °C
Method: ASTM D4052	Measured at 20.0 °C: 778.9 kg/m³
Sample: CycloHEXANE	Corrected at 15.0 °C: 782.8 kg/m³
Commodity group: Refined Product	Status: Test was successful






08-Apr-2025 13:29:44



Density range	0 kg/m³ - 1200 kg/m³
Density resolution	0.01 kg/m³
Density repeatability	0.01 kg/m³ for Distillates and Oils, 0.03 kg/m³ for Gasolines (cold and warm)
Temperature range	10 °C to 40 °C (50 °F to 104 °F)
Temperature stability	0.01 °C (0.02 °F)
Bubble detection method	Automatic Bubble Detection using DensAI
Use	Transparent light and middle distillates

Density CRM Standards:	Volume:	ISO 17034
SETA-1414-0083	250 mL Set of 3	✓

SetaVap2 Automatic Vapour Pressure Tester (81000-2)

ASTM D5191; IP 394; IP 409; IP 481; EN 13016-1

Correlates with ASTM D323; ASTM D4953; D5482; ASTM D6378; IP 69

Pressure range	0 to 200 kPa ±0.5 kPa (0.1 kPa resolution)
Pressure resolution	0.1 kPa
Chamber temperature	37.8 °C (100 °F) 0.1 °C
Modes	Ptot, DVPE, EPA, Crude, Pabs
Power requirement	115 Vac or 220 Vac, 50/60 Hz, 70 W
Size (HxWxD) / Weight	37x13x20 cm / 3.5 kg



SetaCheck® Biodiesel Plus (SA5550-0)


ASTM D8274 (Type 1); IP 639

SetaCheck Biodiesel is a handheld instrument designed to provide fast, accurate and repeatable on-site measurement of biodiesel content in diesel fuel blends.


- Samples measured directly
- Correlates with EN 14078 and ASTM D7371
- Unique Mid-IR technology
- Simple to use
- Bright, backlit screen displays the Biodiesel concentration (% by volume)
- User calibration available via PC software
- 9:1 Dilution mode available for high biodiesel concentrations



Range	0.1% to 80%
Operating temperature	5 - 40 °C
Test duration	Less than 1 minute
Sample volume	2 mL
Power supply	Rechargeable Lithium Ion battery
Size (HxWxD) / Weight	36x15x20 cm / 1 kg

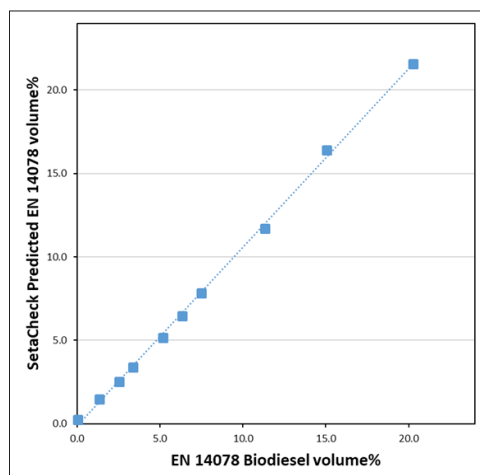

SetaCheck Biodiesel Plus versus EN 14078

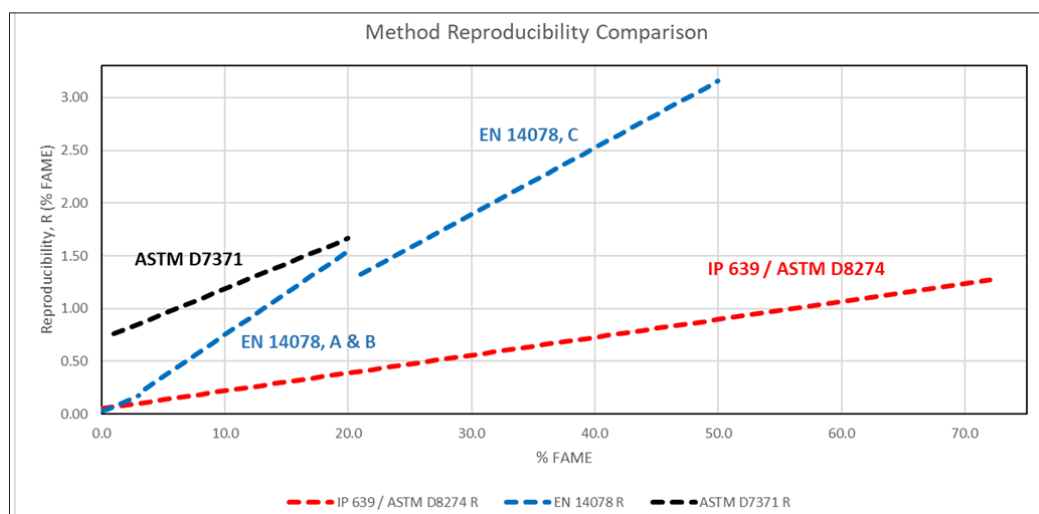

SetaCheck Biodiesel Plus, precision versus other methods



SetaCheck® Ethanol (SA7500-0)

The SetaCheck Ethanol is a portable handheld instrument designed to provide an in-situ measurement of the ethanol content in a gasoline-ethanol fuel blend.

- Portable handheld and easy to use
- 0.1 % to 100 % measurement range
- Less than 10 mL of sample required
- Automatic temperature measurement
- Rapid test time
- Operates on rechargeable battery
- Large colour display with capacitive touchscreen



Range	0.1 to 100 % by volume
Operating temperature	10 to 30°C
Test duration	Less than 1 minute
Minimum sample volume	2 mL
Power	Rechargeable lithium ion battery
Size (H x W x D) / Weight	36 x 15 x 20 cm / 1 kg

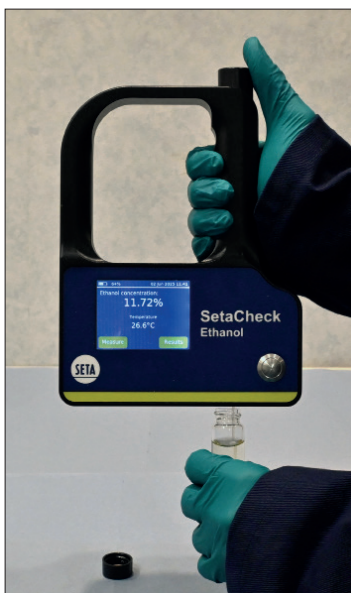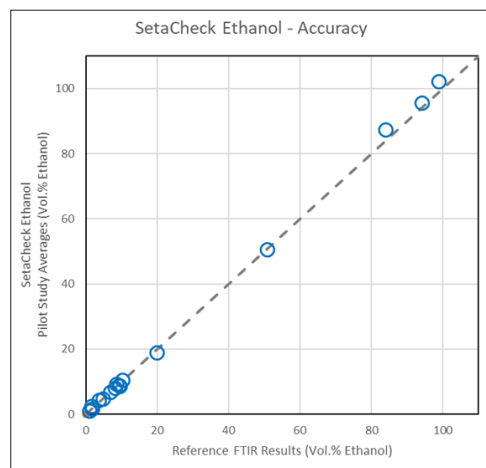

1) Take sample



2) Press measure



3) Test complete

Seta Existent Gum Solid Bath

ASTM D381; IP 131; IP 540; BS EN ISO 6246; DIN 51 784

A Solid Block Bath designed to carry out up to five simultaneous tests for determining existent gum content in fuels by the Jet Evaporation method.

- 5 test stations
- Suitable for air and steam or air only operation
- Integral steam superheater (12200-3 only)
- Digital temperature control
- Flow gauge calibrated for air and steam



Air & Steam Bath (12200-3)

Air Only Bath (12210-0)

	Air & Steam 12200-3	Air only 12210-0
Operating temperature range	140 to 260 °C (± 0.5 °C) (284 to 500 °F)	
Heaters	Integral steam superheater	2 kW
Operational requirement	Steam generator or rotary compressor	Rotary compressor
Over temp cut out	280 °C (adjustable)	
Air/steam inlet	15 mm o.d pipe	
Size (HxWxD) / Weight	45x35x50 cm / 45 kg	



SSAFlab Filtration Kit for ASTM D5452, IP 423 (19722-0)

ASTM D5452; IP 423

SSAFlab Filtration Kit for the determination of particulate contamination in aviation fuels.

Comprises; ground/bond wire, 0.45 micron filter, matched weight pair of 0.8 micron filters, stainless steel filter holder with stainless steel screen, filter forceps, 25 mm solvent filtering dispenser, vacuum filtering flask, petrislide.




Seta Solid Block Bath 4 Way (16670-0)

ASTM D525; ASTM D873; IP 40; BS EN ISO 7536; BS 2000-40; BS 2000-138; IP 138

The Seta Solid Block Bath is an aluminium block bath with digital temperature control that has been designed to accept up to four vessels for the oxidation stability of gasoline and aviation fuels.

- 4 place dry block bath
- Small footprint
- Ambient +5 to 130 °C (heat up time of 80 minutes)
- Stability ± 0.1 °C
- Over temperature cut out
- Auto-Oxi control software available (15452-2)

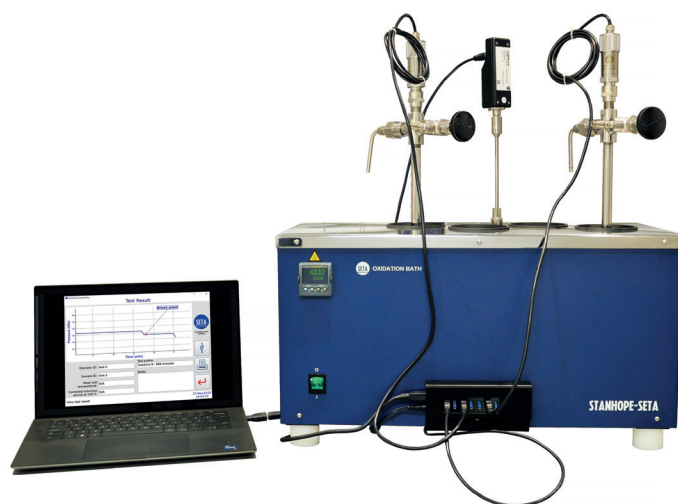

Bath type	Dry
Temperature range	Ambient to 130 °C
Number of test stations	4
Heat up time to 100 °C	80 minutes
Power requirements	200/240 V, 50/60 Hz, 3 kW
Size (HxWxD) / Weight	30x60x28 cm / 65 kg





Pressure vessels include safety burst disc

Seta Auto-Oxi Control System Software (15452-2)

Supports ASTM D525; IP 40; ASTM D873; IP 138; ASTM D5304; ISO 7536

The Seta Auto-Oxi software is a data capture system for simultaneously monitoring up to four pressure vessels.

- Software for use with water or dry block baths
- Incorporates purge and leak test
- Automatic breakpoint determination
- Connects to Seta Pressure Transducers (15455-3)
- Connects to a PC via USB data cable

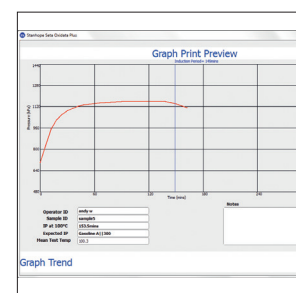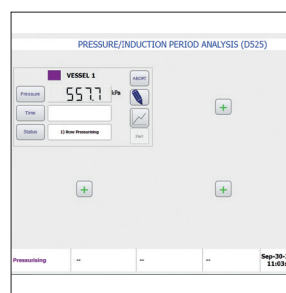

Seta TOST Solid Block Bath (16640-2)

ASTM D943; ASTM D2274; IP 157; IP 388; ISO 4263

The Seta TOST Bath is a compact solid block bath which can accept up to 12 sets of oxidation glassware, it is suitable for TOST tests and as an alternative to liquid filled baths.

- Suitable for TOST tests and as an alternative to liquid filled baths
- 12 positions dry block bath
- Ambient to 100 °C temperature range
- Integral Setatemp controller
- Seta Oxflo controller and flowmeters included



TOST Dry Block Bath - 6 Way with Oxflo (16645-2)

Oxidation Bath Liquid - 6 Way with Oxflo Controller (16900-7)

Temperature range	Ambient to 100 °C
Temperature control at 95 °C	Uniformity ± 0.1 °C Stability ± 0.03 °C
Bath type	Solid block
Power requirements	220/240 V, 50/60 Hz / 2 kW
Size (HxWxD) / Weight	Bath: 50x25x49 cm / 110 kg Oxflo Control Unit: 44x35x26 cm / 11 kg



Seta Rust Prevention Test Bath (11200-8)

ASTM D665; ASTM D3603; ASTM D5534; IP 135; ISO 7120; NACE TM-01-72

The Seta Rust Prevention Test Bath allows determination of the ability of steam turbine oils to prevent rusting of ferrous components and to test the corrosive properties of cargoes in petroleum product pipelines.

- Seven test stations
- Quick release stirrers
- Digital temperature controller
- Bath viewing window and door
- Integral drain valve
- Overflow drain valve



Seta Copper Silver Block Bath (11310-0)

ASTM D130; ASTM D7667

The Seta Copper Silver Corrosion Solid Block Bath features four test wells to accommodate Seta Corrosion Test Vessels (11500-0) and four test tube pre heating wells for 25 mm diameter x 150 mm test tubes.

- Dry heating for up to four test samples
- Accommodates four 25 x 150 mm test tubes and four pressure vessels
- Operator selected test set point temperature
- Stability ± 0.1 °C
- Over temperature cut-off



Seta Silver and Copper Corrosion Bath (11405-0)

ASTM D130; ASTM D4048; ASTM D7667; ASTM D767; IP 112; IP 154; IP 611; BS EN ISO 2160 & BS 2000-154 & ISO 2160; DIN EN ISO 2160; ISO 2160

The Seta Silver and Copper Corrosion Bath is a bench top, stainless steel water or oil bath. Digitally controlled to ± 0.1 °C over a temperature range of ambient +5 °C to 150 °C, the bath is used to determine the corrosive effects of fuels, lubricants and other hydrocarbon products on copper and silver.

- Maximum capacity of 6 test stations
- Precise temperature control
- Digital timer
- Over temperature cut out and liquid level float switch
- Maximum liquid level outlet fitting
- Small footprint



Seta Corrosion Test Bath (22220-5)

ASTM D1838; IP 411; BS EN ISO 6251 & BS 2000-411 & ISO 6251

The Seta Corrosion Test Bath is a stainless steel water bath which accommodates up to 3 pressure vessels. Fitted with heater and thermostatic temperature controller with integral circulator and LED display.

- Temperature range ambient +5 °C to 120 °C
- Stainless steel water bath which accommodates up to 3 pressure vessels
- Fitted with heater and thermostatic temperature controller with integral circulator and LED display



Micro Carbon Residue Tester (97400-3)

ASTM D4530; ASTM D189; IP 13; IP 398; ISO 10370

The Seta MCRT is an automatic instrument designed to determine the carbon residue formed after evaporation and pyrolysis of petroleum products.

- Carbon residue range 0.1% to 30.0% (m/m)
- Fully automatic, load and go
- 12 sample capacity
- Automatic temperature ramp and gas flow control
- Equivalent to ASTM D189; IP 13
- Temperature range ambient to 500 °C
- User friendly interface
- Calibrated for precise temperature ramp rate
- Integrated fan for cooling
- Low pressure and over temperature cut out



Temperature range	Ambient to 500 °C $\pm$ 2 °C
Ramp rate	10 to 15 °C/min
Pressure controller	20 to 750 kPa
Primary pressure	1400 kPa max
Flowmeter	0.1 to 1 litre/min
Flow rate	Automatic, 150 mL/min or 600 mL/min
Size (HxWxD) / Weight	46x35x39 cm / 21 kg



MTVM Fuel Oil Standard:	Test Method:	Range:	Volume:
99852-0	ASTM D4530; IP 398; ISO 10370	0.10 to 30.0 %(m/m)	500 mL

Setaclean Total Sediment Tester (16120-2)

ASTM D4870; BS ISO 10307-1; BS ISO 10307-2;
BS 2000-390; IP 375; IP 390

- Measures sediment up to 0.5% m/m
- Sample viscosities up to 130 m²/s
- 10 g sample size
- Two test stations



SetaCool - Cloud and Pour Point

ASTM D97; ASTM D5853; IP 15; IP 441; BS 2000-15;
ISO 3016; Pour Point of Petroleum Products
ASTM D2500; IP 219; BS EN 23015; BS 2000-219;
ISO 3015; Cloud Point of Petroleum Products
ASTM D6371; IP 309; EN ISO 116; Cold Filter Plugging
Point of Diesel and Heating Fuels

The SetaCool is designed to provide fast, accurate and precise determination of Cloud and Pour Point of petroleum products, diesel and fuels cooling to -78 °C.

- Built in cooling
- Safe for the environment, no F-Gas or CFC's
- Dry block bath with no cooling medium required
- Independent modular system
- Low power consumption and energy loss
- Temperature stability +/- 0.2 °C or better
- Compact footprint
- Low maintenance



*Example configuration

	SetaCool 94155-0	SetaCool-12 94158-0	SetaCool-16 94160-0	SetaCool-20 94165-0
Number of test jackets	4 places x 1 unit	4 places x 3 units	4 places x 4 units	4 places x 4 units
Cooling temperature range	Adjustable between Ambient to -78 °C with PID control			
Temperature stability	+/- 0.2 °C or better			
Power	100 V 240 Vac, 50/60 Hz			
Size (HxWxD) / Weight	Individual unit: 45x32x24.9 cm / 19.4 kg			

Seta Cloud and Pour Point Cryostat

ASTM D97; ASTM D5853; IP 15; IP 441; BS 2000-15;
ISO 3016; Pour Point of Petroleum Products
ASTM D2500; IP 219; BS EN 23015; BS 2000-219;
ISO 3015; Cloud Point of Petroleum Products

- Four individually temperature controlled compartments
- Four air wells in each compartment
- CFC free refrigeration system
- Heated anti-condensation lid

	Seta Cloud and Pour Point Cryostat 93531-9
Temperature range	Ambient to -35 °C (x3) -35 to -51 °C (x1)
Compartment volume	2 Litres
Test positions	16
Cool down time	Approx 2 hrs (from 32 °C to -51 °C)
Power supply	1.8 kW 110/220 Vac, 50/60 Hz
Size (HxWxD) / Weight	95x63x62 cm / 115 kg

Air Release Value System (15840-0)



ASTM D3427; IP 313; BS 2000-313; ISO 9120

Seta Air Release Value (ARV) System is a benchtop, automated instrument, used to evaluate the ability of turbine, hydraulic and gear oils to release entrained air.

- Integrated solution
- Multi-station sample management platform
- Simple test menu
- Automated density monitoring
- Circulated sample heating
- Quick connectors for ease of sample handling
- Integral sinker warmer
- Automatic result calculation
- Results storage for over 10,000 tests
- Full LIMS connectivity



Sample size	200 mL
Test temperature range	Ambient to 75 °C (air to 85 °C)
Set temperatures	25 °C, 50 °C, 75 °C (custom temp available in software)
Sample temperature stability	±0.1 °C
Air temperature stability	±0.2 °C
Water supply	10 l/min, adjustable from 25 °C to 80 °C
Density	0.0001 g/mL (0.1 kg/m ³)
Size (HxWxD) / Weight	82x44x50 cm / 29 kg

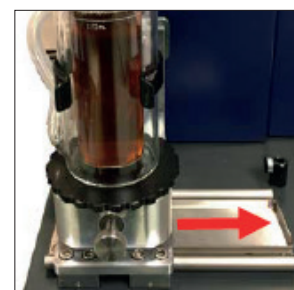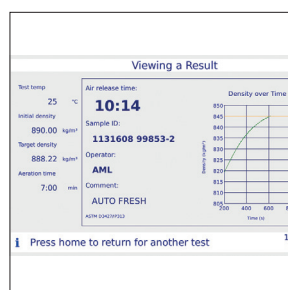

ARV Reference Oil 500 mL (15842-0)

- Produced in general accordance with the principles of ISO 17034
- ASTM D3427 is equivalent to both IP 313 and BS 2000-313
- Certified ARV value at 50 °C and 75 °C
- Valid for 12 months from manufacture
- 14 to 18 minutes at 50 °C
- 5 to 7 minutes at 75 °C



Herschel Emulsifier (96700-2)



ASTM D1401; IP 412; ISO 6614

The Herschel Emulsifier is a bench top, automated instrument, used to measure the ability of petroleum oils or synthetic fluids to separate from water.

- 4 independently controlled test stations
- Multiple operator safety features, including emergency stop button
- Motorised raising and lowering
- Automated test sequence
- Intermediate scraping position
- Integrated timing
- Easy handling and removal of sample
- No removal of the paddles required
- Guaranteed paddle rotation speed
- Non-reflective enhanced LED lighting



Bath volume	5 Litres
Bath liquid	White oil
Sample size	40 mL oil 40 mL distilled water 1% sodium chloride solution or synthetic seawater
Test temperatures	54 °C and 82 °C
Bath temperature stability	±1 °C
Stirrer speed	1500 ±15 rpm
Size (HxWxD) / Weight	89x45x45 cm / 49.5 kg

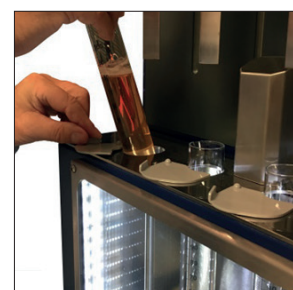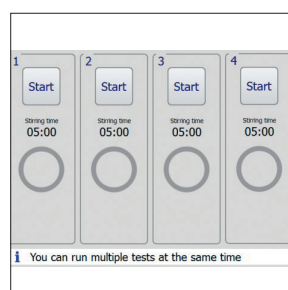

Setafoam Dual Twin Foam Test Baths (14020-8)

ASTM D892; IP 146; ISO 6247

Setafoam Dual Twin Foam Test Baths are a pair of highly transparent water baths for detecting foaming characteristics in lubricating oils.

- High and low temperature baths
- Up to two simultaneous tests per bath
- Two pre-heating stations per bath
- Two integral normalising coils
- Digital or Analogue flow meter options available
- LED back lighting



Spherical Diffuser and Tube (14021-0)

Cylindrical Metal Diffuser and Tube (14023-0)

Stanomatic Dual Grease Workers

ASTM D217; ASTM D7342; BS 2000-50; IP 50; ISO 2137

The Stanomatic Dual Grease Worker removes the variations that can occur in the hand operation of grease workers, particularly with hard greases.

Grease Worker with Brass Cups (17790-3)

Grease Worker with Stainless Steel Cups (17792-3)



Universal Penetrometer (17190-3)

ASTM D5; ASTM D217; ASTM D937; ASTM D1321; ASTM D1403; ASTM D1831; ASTM D7342; IP 49; IP 50; IP 179; IP 310; IP 376

- Penetration range 0 to 640 Pen (64 mm)
- 5, 8, 10, 12, 30 and 60 second timer settings
- Solid state timer
- Automatic plunger release and retention
- Flexi light
- Semi-automatic operation



Dual Water Washout Tester (19650-0)

ASTM D1264; DIN 51 807; IP 215; ISO 11009

The Seta Dual Water Washout Tester estimates the resistance of a lubricating grease to washout by water from a bearing.

A dual bath instrument with the capability for simultaneous testing at 38 °C and 79 °C allowing significant time savings for D1264 tests.

Single Bath Water Washout Tester (19610-4)



Roll Stability Testers

ASTM D1831; ASTM D8022

- 165 rpm or 100 to 200 rpm variable speed control
- Air bath temperature range ambient to 200 °C
- Integral timer
- Variable speed control
- Polished stainless steel or mild steel cylinder/roller options
- Optional venting hood for ambient tests

Seta High Temperature Roll Stability Tester (19400-5)

Seta High Temperature Roll Stability Tester SC, variable speed (19450-0)




KV-10 Viscometer Bath (84300-0)

ASTM D2170; ASTM D2270; ASTM D445; ASTM D446;
IP 71; IP 226; IP 319; BS 188; BS 2000-71; ISO 3105;
BS EN 12595; BS 2000-319; EN 3104; DIN 51 366;
DIN 51 562

The KV-10 Viscometer Bath precisely maintains the temperature of viscometer tubes, which are used to measure the viscosity of liquid petroleum products. The instrument consists of an oil-filled temperature-controlled bath. The bath holds up to ten standard viscometer tubes, which can be viewed through the glass panel at the front of the instrument. Each aperture has a cover to minimise heat loss when not in use. LEDs illuminate the bath so that the viscometer tubes can be viewed more easily.

- Up to 10 viscometer tubes
- 50 litre oil/water bath
- Temperature stability ± 0.02 up to and including 100 °C, ± 0.03 above 100 °C
- Digital display with 0.01 °C resolution
- Oil, silicone fluid or water filled depending on temperature of use
- Low liquid level and over temperature protection
- Double wall glass front panel ensures optimum insulation and reduce heat loss
- Integral back lighting
- Clear and precise reading of viscometer tubes

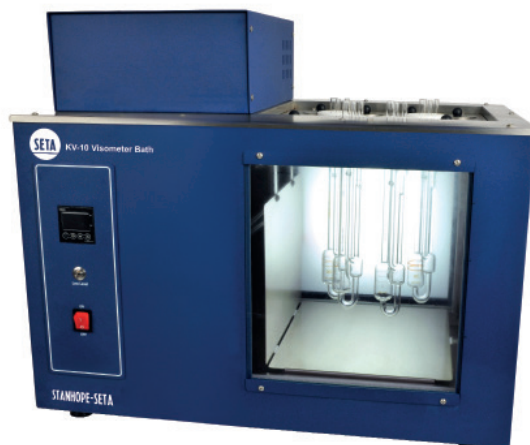

Temperature range	Ambient to 150 °C (302 °F)
Temperature stability	± 0.02 °C up to and including 100 °C, ± 0.03 °C above 100 °C
Tube capacity	up to 10
Bath fluid	Oil/silicone/water
Bath capacity	40 Litres
Power requirements	220/240 V, 50 Hz / 2.2 kW
Size (HxWxD) / Weight	41x63x38 cm / 25 kg

KV-2 Viscometer Bath (94710-4)

ASTM D445; ASTM D446; IP 71; BS 188; BS EN ISO 3104;
BS 2000-71.1; ISO 3104; BS 2000-71-2; ISO 3105

- Up to 2 viscometer tubes
- Temperature range -40 °C to +20 °C
- High stability ± 0.01 at -20 °C
- Illuminated bath
- Insulated and heated viewing window
- Heated top plate

Temperature range	-40 °C to 20 °C
Temperature stability	± 0.01 °C at -20 °C
Tube capacity	up to 2
Bath fluid	Anhydrous methanol or water/ ethylene glycol (50/50)
Bath capacity	7 Litres
Size (HxWxD) / Weight	75x37x70 cm / 56 kg



Salt in Crude Analyser (99700-6)

ASTM D3230; IP 265

The Seta Salt in Crude Analyser is a robust and portable instrument for determining the chloride (salt) content of crude oils.

- Automatic salt concentration calculation
- Pre-calibrated for immediate use
- No need for mixing of salt standards
- Typical test time of less than 30 seconds
- Interchangeable sensor
- 3 operating modes (ASTM, IP, user defined)
- Results can be displayed or exported to PC or a LIMS network
- Battery or mains voltage operation
- Moistureproof
- Fully portable



Conductivity range	0.0 to 151 lb / 1000 bbl (Res: 0.1 lb/bbl) 0.0 to 430.0 g/m ³ (Res: 0.1 g/m ³)
Temperature range	-20 °C to 150 °C (Res: 0.1 °C)
Power supply	9 Vdc battery or mains adaptor 100/220 Vac, 50/60 Hz
Size (HxWxD) / Weight	23x12x5 cm / 1.6 kg



SaltCheck Verification Tool (99703-0)

ASTM D3230; IP 265

A set of 3 'SaltCheck' verification modules, supplied in a convenient storage case. The modules are calibrated to show equivalent salt values of:

- 0 g/m³ | 30 g/m³ | 190 g/m³



ASTM D3230, Mixed Salt Solution, 250 mL (99704-001)

ASTM D3230, Neutral Refined Oil, 250 mL (99704-002)



Seta Oil Test Centrifuge for Sediment and Water

ASTM D91; D893; D1290; D1796; D2273; D2709; D2711; D4007; ISO 3734; ISO 9030; ISO 15301

The Seta Oil Test Centrifuge is a heating centrifuge fitted with a rotor head with four or six universal pivoting buckets, used to determine water and sediment in oils.

- Heated chamber, ambient to 80 °C
- Static, near vertical bucket positioning
- 4 or 6 place swing out rotor
- Pre-heat facility to ensure bowl is at test temperature
- Conforms to IEC 1010-1 & 1010-2-D
- Optional vapour extraction system



4-Place Centrifuge (90000-3)

6-Place Centrifuge (90100-0)

Maximum RCF	2214 g
Maximum RPM	3000
Temperature range	20 to 80 °C (± 1 °C) (86 to 176 °F)
Power requirement	230 Vac (110 V transformer option) 50/60 Hz / 1 kW
Size (HxWxD) / Weight	41x56x64 cm / 72.5 kg



Crude Oil Filtration Apparatus (19727-0)

ASTM D4807

The Crude Oil Filtration Apparatus determines the sediment in crude oils by membrane filtration in accordance with ASTM D4807. This test method has been validated for crude oils with sediments up to approximately 0.15 mass %.



Filtration Kit for Middle Distillates, IP 440 & EN 12622 (19720-0)/ASTM D6217 (19724-0)

Filtration Kit for Biodiesel & Biodiesel Blends, ASTM D7321 (19730-0)

Seta JetDC (88500-0)

IP 638; ASTM D924; D4052; ASTM D4054; IP 365; D7826

The Seta JetDC measures the dielectric constant and temperature relationship of aviation fuel to predict behaviour in aircraft gauging systems.



- Automated temperature scanning method
- Rapid test (30 minutes per temperature)
- Integrated apparatus with density
- Automated measurement
- Large touchscreen
- Low temperature option with chiller

Sustainable Aviation Fuel (SAF) represents an important route to reduce global net CO₂ emissions. SAF is produced through one of many different processes, using different feedstock such as used cooking oil, syngas, fats, vegetable oils, greases, sugars and alcohols as specified in ASTM D7566. SAFs are required to be approved for use through a standardised testing process, following ASTM D4054, which requires dielectric constant.

Due to a subtle difference in the molecular make up of a traditional Jet Fuel and SAF, there can be a difference in the dielectric constant. Dielectric constant is used by aircraft gauging systems as part of the mechanism to determine how much fuel is on board the aircraft. Aircraft fuel gauge accuracy relies on the dielectric constant properties to be similar across jet fuel batches.




Dielectric constant (κ) measurement range	1.0 to 2.5
Density range (kg/m ³)	1 to 900
Control temperature	0 °C to 40 °C
Power	110/240 V, 50/60 Hz, 300 W
Size (HxWxD) / Weight	40x26x34 cm / 15 kg

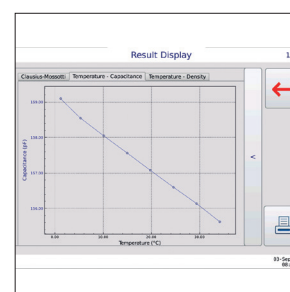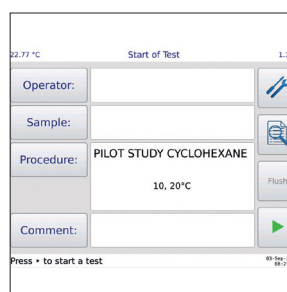

Standardisation

An Energy Institute Task Group has developed a new standard test, IP 628, specifically for aviation fuel which addresses the technology and needs gap of the existing method, ASTM D924 which was identified in CRC report AV-29-20.

IP 638 Determination of Relative Permittivity (Dielectric Constant) of Aviation Fuel, Small Scale Automated Temperature Scanning Method was approved and published in 2023 with full precision. It was developed from proposed method IP PM-FC.

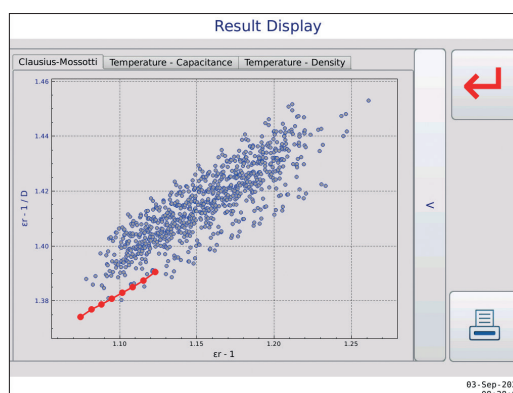

FIJI FAME Measurement

ASTM D7797; ASTM D7963; IP 583;
Specification compliant for ASTM D1655; Defence Standard 91-091; CAN/CGSB-3.23; AFQRJOS; MIL-DTL-83133;
ISO 8217; ISO/PAS 23263

The FIJI FAME Analysers offers full compliance with industry approved methods for rapid and easy checks on parts per million levels of fatty acid methyl ester (FAME) content in fuels, C8 to C22.

- Analysis time 20 minutes
- Certified range 10-150 mg/kg FAME in AVTUR
- Measuring range 0-1000 mg/kg in AVTUR
- Certified range 20 mg/kg to 50% in Distillate and Residual fuels
- Approved for co-processing
- Fully automatic
- Suitable for untrained operators
- No cleaning solvents required
- No pre-sample preparation required
- 50 mL sample volume

FIJI FAME in Jet (SA5000-2)

FIJI FAME in Middle Distillate & Residual Fuel (SA5100-0)


FIJI FAME MultiFuel (SA5200-0)

Sample size	50 mL
Connectivity	2 x USB
User interface	Colour touch screen
Power supply	50 W
Size (HxWxD) / Weight	55x38x42 cm / 27 kg

FIJI uniquely uses state of the art FTIR (Fourier Transform Infra-Red Spectroscopy) technology coupled with a patented sample preparation system. FAME detection is from 10 mg/kg and the instrument has full ASTM and IP method precision from 10 mg/kg level upwards.

FIJI can be used as a lab or field based screening tool. Rapid screening of fuels using the FIJI technique avoids the need for both expensive and complex analytical tests, and can prevent costly delays when releasing jet fuel.

The FIJI instrument is robust, extremely easy to use and is fully automatic so no specialist operator training is involved. Tests require less than 50 mL of sample and typically take under 20 minutes. Results are presented in mg/kg units together with an optional traffic light system for flagging go/no go FAME contamination levels of the fuel.

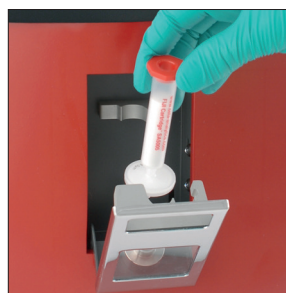


1) Load cartridge



2) Place sample



3) Press GO

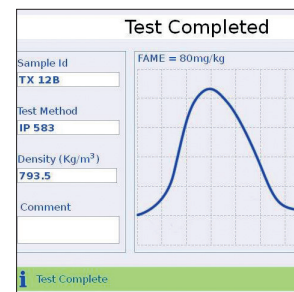

4) Test completed

FAME CRM Standards:	Nominal Value:	Volume:	ISO 17034
SETA-0112-0004	30 mg/kg	100 mL	✓
SETA-0112-0005	100 mg/kg	100 mL	✓

AFIDA (SA6000-0)

ASTM D8183; IP 617; EN 17155

ASTM D4054; ASTM D975; ASTM D6751; ASTM D7467; EN 590; EN 15940; ISO 8217
Correlates with ASTM D613; EN ISO 5165 & DIN 51773

AFIDA provides a fully automated determination of the Indicated Cetane Number (ICN) of diesel and diesel related fuels.

- Indicated Cetane Number by CVCC*
- Highly automated with auto sampler
- Reference Fuel Calibration
- No bias correction to engine values
- Excellent precision to minimise giveaway
- Easy to use, suitable for 24/7 operation
- Load and go

Industry
Approved
Specification
Compliant



*Constant Volume Combustion Chamber

Carousel capacity	36 samples
Range for ICN	35 - 85
Chamber temperature	580 °C
Chamber pressure	17.5 bar
Injection pressure	1000 bar
Compressed air	20.9 ± 0.5 % O ₂
Sample volume	Approx 40 ml for analysis and cleaning
Warm up time	Approx 45 mins
Analysis time	Approx 25 min per sample
Size (HxWxD) / Weight	80x130x60 cm / 160 kg (100 kg excluding platform)

AFIDA incorporates a unique and patented high pressure injection system that generates fine fuel droplets similar to modern common rail injectors in most diesel engines.


A temperature controlled piezo electric injector provides highly repeatable fuel metering, offering improved performance and consistency when compared with solenoid controlled injectors and pintle type nozzles.

The analyser provides very fast, efficient and calibrated ICN determinations using an integral 36 position carousel and auto sampler.


AFIDA Research unit is also available with options for variable chamber temperature, chamber pressure and injection pressure. Gasoline injection system can be installed to give an indication of Octane Number values. For more information, please contact Stanhope-Seta.




1) Place sample in carousel and select location

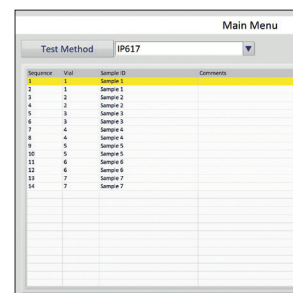

2) Select test sequence and press go

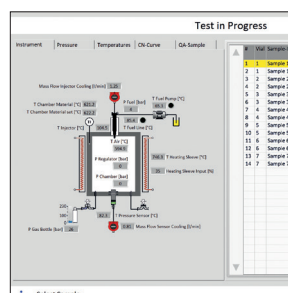

3) Monitor instrument status, pressure, temperature or CN curve during test

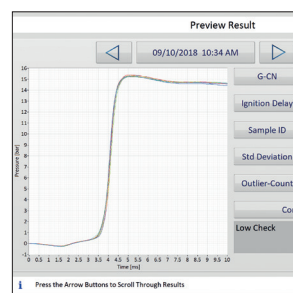

4) View results

H₂S Analyser with Vapour Phase Processor (SA4000-4)

ASTM D7621; IP 570 Procedure A
ISO 8217; ISO/PAS 23263



The H₂S Analyser is a compact bench-top instrument, used to measure the total hydrogen sulphide (H₂S) content of fuel oils and dead crude oils.

- Measurement range from 0-250 mg/kg H₂S (0-250 ppm H₂S) in the liquid phase
- Built in Vapour Phase Processor proven to eliminate effects of chemical interference
- IP 570 Procedure A, referee test method in ISO 8217
- Reference test method IP 570, Procedure A
- A non chemical method, no wet chemistry involved
- Small lab bench footprint
- Integral camera to aid sample input
- Custom result reporting and LIMS output
- Automated gas calibration
- Suitable for monitoring residual marine fuel blends, refinery feedstock components, cargoes and products in the distribution system as well as crude oils
- Critical measurement method for product safety and release

Industry
Approved
Specification
Compliant



Measurement range	0-250 mg/kg
Viscosity range	1-3000 mm ² /s
Principle of measurement	Electro-chemical sensor
Test duration	15 minutes
Sample size	1 mL, 2 mL, 5 mL (depending on H ₂ S concentration)
Diluent volume	20 mL
H ₂ S size (HxWxD) / Weight	52x29x57 cm / 20 kg

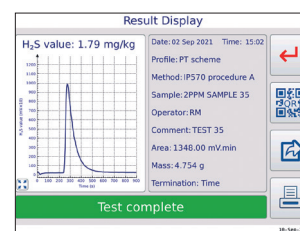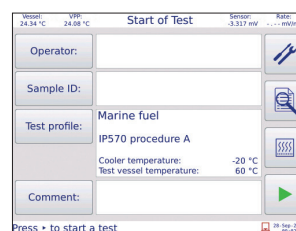

IP 570 Proficiency Test Scheme (PTS) (SA4032-0)

ASTM D7621; IP 570

- Evaluate and monitor laboratory and instrument performance
- Demonstrate compliance with laboratory accreditation requirements
- World-wide inter-laboratory comparison of test results and performance
- Identify any potential equipment or operational bias
- Provide added confidence to laboratory staff and customers



Seta AvCount3 (SA1100-0)

Fuels • Hydraulic oil • Light lubricants

IP 565; IP 630; ASTM D7619; ASTM D6786; ASTM D7647; ISO 11500;
ISO 4406; NAS 1638; SAE AS 4059; GOST 17216-71; NAV AIR 10-1A-17; GJB 420A;
ISO 11171; IP PM-FI
ASTM D975; Defence Standard 91-091; Defence Standard 91-86; AFQRJOS; EN 590; ASTM D1655; ASTM D7467

The AvCount3 is a compact bench-top automatic particle counter, used to measure the size and distribution of particles and water droplets in light and middle distillate fuels, including aviation fuel and kerosine, biodiesel, low viscosity oils and hydraulic oils.

- ISO 11171 calibration
- Cumulative counts/ml
- ISO 4406 Cleanliness Codes
- Colour touch screen
- Dilution ratio calculation
- Real time display of test progress
- User programmable
- Under 3 minute test time (IP 565)
- LIMS, network and VNC connectivity
- Programmable alarm limits
- User and sample identification
- 14 embedded test methods
- Integrated printer
- 500,000 test memory
- In-field verification and calibration



Start of Test

Operator:		
Sample:		
Method:	ISO 4406	
Dilution:	0%	No multiplier
Comment:		

1 measurement Sample: 200 ml Diluent: None

Press ▶ to start a test

- 1) Enter operator and sample details, select method, press ▶

Test in Progress

Flushing

Flush 10ml

Measure 1

4µm cumulative counts per ml (corrected): ---

Sample: FUEL A Operator: RAA Method: ISO 4406 Dilution: 0%

- 2) Test begins, instrument sequences are detailed

Result Display

21 Jun 2021 15:41

Operator: RAA

Sample: FUEL A

Method: ISO 4406

Dilution: 0% (No multiplier)

Comment:

Corrected	Raw data
Size (µm)	Particles per ml
4.0	8298.6
6.0	1144.9
14.0	65.8
25.0	10.4
38.0	5.6
70.3	3.4

Test passed 20/17/13

- 3) Final result displays either numerically or graphically

Counts per measurement (max)	60,000
Sample viscosity (max)	68 mm ² /s
Test methods	14 embedded, user programmable
Number of measuring channels	16
Results	500,000 measurements internal, unlimited on PC Print via internal printer, export to LIMS, USB or QR code
Sample temperature range	0 to 70 °C
Connectivity	RJ45 Ethernet or USB
Power supply	100 – 240 Vac, 50/60 Hz auto-sensing universal power supply
Size (HxWxD) / Weight	37x23x27 cm / 6 kg

Seta AvCountHV (SA1110-0)

Hydraulic oil • lubricants

ASTM D6786; ASTM D7647; ISO 4406; ISO 11171; ISO 60970; NAS 1638

The AvCountHV is a compact bench-top automatic particle counter, used to measure the size and distribution of particles and water droplets in lubricant and hydraulic oils with viscosities up to 1500 mm²/s at 20 °C.

- High viscosity samples
- ISO 11171 calibration
- Cumulative counts/ml
- ISO 4406 Cleanliness Codes
- Colour touch screen
- Dilution ratio calculation
- Real time display of test progress
- User programmable
- LIMS, network and VNC connectivity
- Programmable alarm limits
- User and sample identification



Featuring a simple user interface with touch screen display, the fully-automated test means extensive operator training is not required before using the instrument. The simplicity of test set up and automation, gives operators the time to work on something else, helping to reduce labour costs.

The AvCountHV's enhanced functionality includes a dilution protocol for testing of high viscosity samples, real-time display of test progress, 4 embedded test methods, user and sample identification and password protected levels.

Statistical Quality Control software allows analysis of results in accordance with ASTM D6299.

Counts per measurement (max)	30,000
Sample viscosity (max)	1500 mm ² /s
Test methods	4 embedded, user programmable
Number of measuring channels	16
Results	500,000 measurements internal, unlimited on PC Print via internal printer, export to LIMS, , USB or QR code
Sample temperature range	0 to 70 °C
Connectivity	RJ45 Ethernet or USB
Power supply	100 – 240 Vac, 50/60 Hz auto-sensing universal power supply
Size (HxWxD) / Weight	37x23x27 cm / 6 kg

SETA-CRM AvCount Verification Material 250 mL (SA1006-0)

- ISO 17034 UKAS accredited
- Manufactured in accordance with ISO 11171 Annex F
- Supplied in 250 mL bottles with a Certificate of Measurement
- Certified for ASTM D7619; IP 565; IP 630



Seta AvCount Lite (SA1800-2)

Fuels • Hydraulic oil • Light lubricants

Defence Standard 91-091; AFQRJOS; EN 590; ASTM D1655; ASTM D975; IP 565; IP 630; ASTM D7619; ASTM D6786; ASTM D7647; ISO 11500; ISO 4406; NAS 1638; SAE A5 4059; GOST 17216; NAVAIR; GJB 420A; ISO 11171; IP PM-FI

The AvCount Lite particle counter provides reliable results for determining the particle concentration in liquid fuels and oils.

Industry
Approved
Specification
Compliant



- ISO 11171 calibration
- Cumulative Particles/ml
- ISO 4406 Cleanliness Codes
- Simple operation
- Portable, compact instrument
- Under 3 minute test time (IP 565)
- Bottle samples or on-line (high and low pressure)
- Stand alone or PC controlled (ProTrend SA1810-0)
- Programmable via PC
- Integrated printer
- Battery power optional

Counts per measurement (max)	60,000
Sample viscosity (max)	64 mm ² /s
Test methods	3 stored at any one time, unlimited on PC
Number of measuring channels	6 displayed (16 recorded and displayed via PC)
Results	900 measurements internal, unlimited on PC
Sample temperature range	0 to 70 °C
Connectivity	USB connection for PC
Power supply	100 – 240 Vac, 50/60 Hz power adaptor Optional 12 Vdc battery and charger
Size (HxWxD) / Weight	25x14.5x35 cm / 7.5 kg



> Online sampling



AVCOUNT LITE	
21/05/18	D7619
16:05:59	IP565
	I4406
	Results
	Service
M: 24/600	

4.0u	248.5	15	D7647
6.0u	68.5	13	
10u	19.0	11	
14u	8.3	10	
21u	1.8	8	
38u	0.3	5	Menu
Meas: Result			

Particle Counting CRM Standards:	Volume:	ISO 17034
SA1006-0	250 mL	✓


Fuels • In-service lubricants
Hydraulic fluids • Insulating oils

The AvCount Lube is a fully configured particle counter and sample delivery system suitable for testing higher viscosity samples such as lubricating oils.


- ISO 11171 calibration
- Up to 200 mm²/s viscosity
- FFKM seals resistant to many synthetic oils
- PC controlled via ProTrend software
- User friendly software
- Integral compressor
- Programmable via PC
- Under 3 minute test time



Counts per measurement (max)	60,000
Sample viscosity (max)	400 mm ² /s with integral Sample Delivery System
Total sample volume used	80 mL (IP 565)
Number of measuring channels	6 displayed (16 recorded and displayed via PC)
Results	600 measurements internal, unlimited on PC
Sample temperature range	0 to 70 °C
Connectivity	USB connection for PC
Power supply	100 – 240 Vac, 50/60 Hz
Size (HxWxD) / Weight	64x32x28 cm / 16 kg

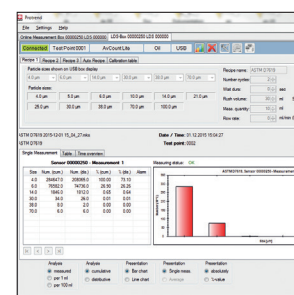

The test process is fully automated. Once a sample is loaded, a test method selected and the test is initiated, the test proceeds without any further operator intervention.

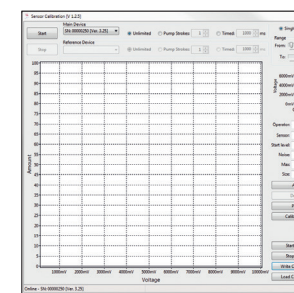[illegible]

Test results are analysed and reported using the ProTrend software, in accordance with a number of test methods and cleanliness codes. For NAS 1638 tests, the sizes are converted automatically using industry accepted conversion tables.

Particle Counting CRM Standards:	Volume:	ISO 17034
SA1006-0	250 mL	✓

Ultrasonic Bath to ISO 11171 (99320-4)

ISO 11171


- Compact benchtop design
- Quiet and safe operation
- Robust and reliable



AvCount 250ml Bottle Tumbler (SA1010-0)

ASTM D975; ASTM D7619; ASTM D7647; IP 565; AS 4059;
Defence Standard 91-091; Defence Standard 91-86; GB 5930GJB

- Fitted with 250 ml 'syrup' bottle holder (SA1004-0)

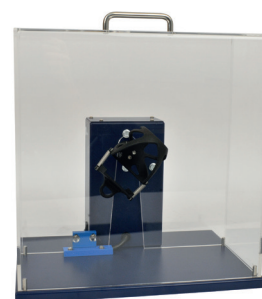

AvCount Calibration Material 'd' 250 mL (SA1121-0)

- ISO 17034 UKAS accredited
- Traceable to NIST 2806d
- For use with particle counters calibrated in accordance with ISO 11171
- Supplied in 250 mL bottles with a Certificate of Measurement



SETA-CRM AvCount Verification Material 250 mL (SA1006-0)

- ISO 17034 UKAS accredited
- Manufactured in accordance with ISO 11171 Annex F
- Supplied in 250 mL bottles with a Certificate of Measurement
- Certified for ASTM D7619; IP 565; IP 630



Water Separation Instrument (WSI) (SA9000-0)

ASTM D8073; IP 624

ASTM D1655; JIG Bulletin No 150; AFQRJOS; ATA 103; CAN/CGSB-3.23 2020; Defence Standard 91-091

The WSI measures how effectively a fuel sample releases entrained and emulsified water when pumped through a water coalescing filter.

- Substantial cost saving per test
- Fully automated sample handling
- Simple operation for non-chemists
- Patented Ultra Sonic Mixer
- 10 minute test duration
- Minimal operator time required (under 3 minutes)
- Sonic mixing sequences providing consistent water droplet size and stable emulsion
- Featuring an API/EI 1581 5th edition filter material
- High resolution optical water detector measuring to <1ppm water concentration



Measurement range	5-100 Water Separation Index (WSI)
Temperature range	18 °C to 29 °C
Sample size	230 mL
Power	Universal, AC 85-264 Vac, 50-60 Hz
Outputs	USB, digital display
Size (HxWxD) / Weight	38x31x35 cm / 8 kg

The WSI is designed to predict filter coalescer failure. The instrument is fully automatic with simple touch screen operation. The WSI features a patented Sonic Mixing Technology which eliminates any variables with water droplet size and ensures the same amount of energy is applied to each sample.

Filter cartridges use latest generation API/EI 5th edition filter material which is less sensitive to weak surfactants and eliminates uncertainty over static dissipator additives, but is responsive to the presence of ester based additives.

In use, an automated test sequence first emulsifies the sample with dyed water using the ultra-sonic mixer. Once emulsified the sample is pumped through the filter cartridge and onto the sensitive water detection system. As the test progresses the filter effectiveness decreases and the detection system measures the water passing through the filter. At the end of the test the result is shown as the Water Separation Index (WSI).



- WSI Replacement Filters (pack of 50) SA9001-008
- Used in conjunction with the Filter Assembly Tool SA9015-0
- 2 year shelf life from date of manufacture



- WSI Pre-assembled Filter (pack of 100) SA9001-006
- Load and go
- 2 year shelf life from date of manufacture

WSI RM Standards:	Range:	Volume:
SA9004-0	97.5 to 100 Water Separation Index	100 mL

Multi Test Verification Materials

- Highly cost effective solution to laboratory verification requirements
- Verification to ASTM/CEN/ISO/IP
- Certified values
- 2 year shelf life from manufacture

99850-0 Seta MTVM - Kerosene (Jet Turbine Fuel) 500 mL

Test Name	ASTM / Method	Range	Amount/ test
Density 15 °C	D4052; IP 365; ISO 12185	0.798 to 0.807 kg/l	2 mL
Distillation IBP	D86; IP 123; EN ISO 3405	140 to 180 °C	100 mL
Distillation 10%	D86; IP 123; EN ISO 3405	159 to 188 °C	100 mL
Distillation 50%	D86; IP 123; EN ISO 3405	183 to 218 °C	100 mL
Distillation 90%	D86; IP 123; EN ISO 3405	220 to 252 °C	100 mL
Distillation FBP	D86; IP 123; EN ISO 3405	241 to 269 °C	100 mL
Distillation Residue	D86; IP 123; EN ISO 3405	1.1 to 1.3%(v/v)	100 mL
Distillation Loss	D86; IP 123; EN ISO 3405	0.4 to 0.8%(v/v)	100 mL
Flash Point	IP 170; EN ISO 13736	35 to 60 °C	85 mL
Freezing Point	D2386; D7153; D7154; IP 16	-62 to -43 °C	25 mL
Aromatics HPLC	D6379; IP 436	0.7 to 50%(v/v)	2.5 mL
Smoke Point	D1322; IP 57; ISO 3014	20 to 25 mm	20 mL
Kin Vis -20 °C	D445; IP 71; EN ISO 3104	2 to 8 mm ² /s	20 mL
Acid Number	D3242; IP 354	< 0.100 mg KOH/g	100 mL
Mercaptans	D3227; IP 365; ISO 12185	0.0003 to 0.0100% (m/m)	40 mL

99852-0 Seta MTVM - Fuel Oil 500 mL

Test Name	ASTM / Method	Range	Amount/ test
Density at 15 °C	D4052; IP 356; ISO 12185	0.94 to 0.994 kg/L	2 mL
Pour Point	D97; IP 15; EN ISO 3016	-14 to 17 °C	Up to 45 mL
Kin Vis 50 °C	D445; IP 71; EN ISO 3104	150 to 1800 mm ² /s	Up to 40 mL
Kin Vis 100 °C	D445; IP 71; EN ISO 3104	20 to 95 mm ² /s	Up to 40 mL
Micro Carbon	D4530; IP 398; ISO 10370	0.10 to 30.0%(m/m)	2 mL
Flash Point	D93 (b); IP 34 (b); EN ISO 2719 (b)	92.3 to 121.6 °C	75 mL

99856-0 Seta MTVM - Bitumen 200mL

Test Name	ASTM/Method	Range	Amount/ test
Softening Point	IP 58; EN ISO 1427	37 to 54 °C	7.5 mL

99851-0 Seta MTVM - Gas Oil 500 mL

Test Name	ASTM / Method	Range	Amount/ test
Density at 15 °C	D4052; IP 356; ISO 12185	0.83 to 0.854 kg/L	2 mL
Distillation IBP	D86; IP 123; EN ISO 3405	160 to 190 °C	100 mL
Distillation 10%	D86; IP 123; EN ISO 3405	198 to 242 °C	100 mL
Distillation 50%	D86; IP 123; EN ISO 3405	255 to 290 °C	100 mL
Distillation 90%	D86; IP 123; EN ISO 3405	319 to 350 °C	100 mL
Distillation 95%	D86; IP 123; EN ISO 3405	335 to 368 °C	100 mL
Distillation FBP	D86; IP 123; EN ISO 3405	347 to 385 °C	100 mL
Distillation Residue	D86; IP 123; EN ISO 3405	1.25 to 1.42%(v/v)	100 mL
Distillation Loss	D86; IP 123; EN ISO 3405	0.26 to 0.7%(v/v)	100 m:
FAME Content	IP 579	0.03 to 19.03%(v/v)	
Flash Point	D93; IP 34; EN ISO 2719	56 to 80 °C	75 mL
Cloud Point	D2500; D5771- IP 444; D5772; IP 445; D5773; IP 446; IP 219; ISO 3015; EN 23015	-17 to -4 °C	Up to 38 mL
CFPP	D6371; IP 309; EN 116	-30 to 0 °C	45 mL
Pour Point	D97; IP 15; D5950; D5949; D6749; D6892; D5985; EN ISO 3016	-33 to -6 °C	Up to 38 mL
Kin Vis 40 °C	D445; IP 71; EN ISO 3104	2.3 to 3.5 mm ² /s	Up tp 40 mL
Lubricity HFRR	D6079; IP 450; ISO 12156-1	212 to 512 µm	2 mL
Water Karl Fischer	D1744; IP 438; EN ISO 12937	23.4 to 63.9 mg/kg	5 mL

99854-0 Seta MTVM - Motor Gasoline 500 mL

Test Name	ASTM / Method	Range	Amount/ test
Motor Octane No.	D2700; ISO 5163; IP 236	82 to 90 MON	-
Research Octane No.	D2699; ISO 5164; IP 237	90.5 to 101 RON	-
Density at 15 °C	D4052; IP 356; ISO 12185	0.763 kg/L	2 mL
Distillation IBP	D86; IP 123; EN ISO 3405	26 to 39 °C	100 mL
Distillation 70°C	D86; IP 123; EN ISO 3405	5 to 44%(v/v)	100 mL
Distillation 100°C	D86; IP 123; EN ISO 3405	36 to 70%(v/v)	100 mL
Distillation 150°C	D86; IP 123; EN ISO 3405	75 to 95%(v/v)	100 mL
Distillation FBP	D86; IP 123; EN ISO 3405	175 to 205 °C	100 mL
Motor Gasoline Aromatics	D6839; ISO 22854; EN 14517	19.32 to 46.29%(v/v)	0.75 mL
Motor Gasoline Olefins	D6839; ISO 22854; EN 14517	0.40 to 26.85%(v/v)	0.75 mL
Motor Gasoline Saturates	D6839; ISO 22854; EN 14517	26.85 to 79.31%(v/v)	0.75 mL
Vapour Pressure	D5191; IP 394; EN 13016-1	50 to 95 kPa	1-10 mL

Single Test Verification Materials

- Verification to ASTM/CEN/ISO/IP
- Certified values
- 2 year shelf life from manufacture

Seta STVM - Kerosene 50 mL

Seta Part Number	Test Name	ASTM/Method	Range
99857-0	Napthalenes	D1840	0.08 to 5.6% vol
99858-0	Sulfur ED X-ray	D4294; IP 336; EN ISO 8754	0.003 to 0.07%(m/m)
99898-0	Aromatics-Di HPLC	D6379; IP 436	1.44 to 2.82%(m/m)
99898-0	Aromatics-Mono HPLC	D6379; IP 436	15.2 to 21.2%(m/m)
99898-0	Aromatics-Total HPLC	D6379; IP 436	16.9 to 22.8% m/m

Seta STVM - Gas Oil 50 mL

Seta Part Number	Test Name	ASTM/Method	Range
99867-0	Sulfur WD X-ray	D2622; ISO 20884	0.001 to 0.013%(m/m)
99869-0	Aromatics HPLC	D6591; IP 391; EN 12916	0.6 to 3.6%(m/m)

99905-0 Seta SIMDIS STVM - Kerosene (Jet Turbine Fuel) 10 mL

Test Name	ASTM/Method	Range	Amount/Test
SIMDIS IBP	D2887; IP 406	100 to 180 °C	1 mL
SIMDIS 10%	D2887; IP 406	142 to 188 °C	1 mL
SIMDIS 50%	D2887; IP 406	183 to 218 °C	1 mL
SIMDIS 90%	D2887; IP 406	220 to 249 °C	1 mL
SIMDIS 95%	D2887; IP 406	230 to 259 °C	1 mL
SIMDIS FBP	D2887; IP 406	244 to 297 °C	1 mL

99906-0 Seta SIMDIS STVM - Gas Oil 10 mL

Test Name	ASTM/Method	Range	Amount/Test
SIMDIS IBP	D2887; IP 406	120 to 190 °C	1 mL
SIMDIS 10%	D2887; IP 406	177 to 242 °C	1 mL
SIMDIS 50%	D2887; IP 406	260 to 290 °C	1 mL
SIMDIS 90%	D2887; IP 406	320 to 355 °C	1 mL
SIMDIS 95%	D2887; IP 406	335 to 375 °C	1 mL
SIMDIS FBP	D2887; IP 406	350 to 414 °C	1 mL

99980-0 Seta STVM - Gasoline 5 x 5 mL

Test Name	ASTM/Method	Range	Amount/Test
Motor Gasoline Aromatics	D6839; ISO 22854; EN 14517	22.41 to 35.41%(v/v)	0.1 µL
Motor Gasoline Olefins	D6839; ISO 22854; EN 14517	0.31 to 14.50%(v/v)	0.1 µL
Motor Gasoline Saturates	D6839; ISO 22854; EN 14517	43.62 to 77.2%(v/v)	0.1 µL
Motor Gasoline Oxygenates	D6839; ISO 22854; EN 14517	0.04 to 12.3%(v/v)	0.1 µL
Motor Gasoline Oxygen	D6839; ISO 22854; EN 14517	0 to 2.58%(v/v)	0.1 µL
Motor Gasoline Benzene	D6839; ISO 22854; EN 14517	0.051 to 0.888%(v/v)	0.1 µL

99981-0 Seta STVM - Gas Oil 50 mL

Test Name	ASTM/Method	Range	Amount/Test
Gas Oil Sulphur	D5453; IP 490; ISO 20846	3 to 10 mg/kg	1 mL

Worldwide service and support is available through Stanhope-Seta's network of experienced and trained representatives.

To find your nearest representative, visit www.stanhope-seta.co.uk/local-representative/



Contact Us

Stanhope-Seta
London Street
Chertsey
Surrey
KT16 8AP
UK

t: +44 (0) 1932 564391
e: sales@stanhope-seta.co.uk

www.stanhope-seta.co.uk
www.seta-crm.co.uk

